


Repositionable Breast Localization Needles

The Homer Repositional Breast Localization Needle is a strong, flexible J-curve wire that may be retracted and repositioned as often as necessary to achieve accurate placement.

- J-curve wire and stabilizer design locks the needle and wire in place, anchoring the system to prevent wire migration while protecting against unguarded needle penetration during breast compression and patient transport
- System provides a palpable guide to the surgeon, thereby offering a choice of paths to the lesion
- Centimeter markings facilitate accurate depth placement

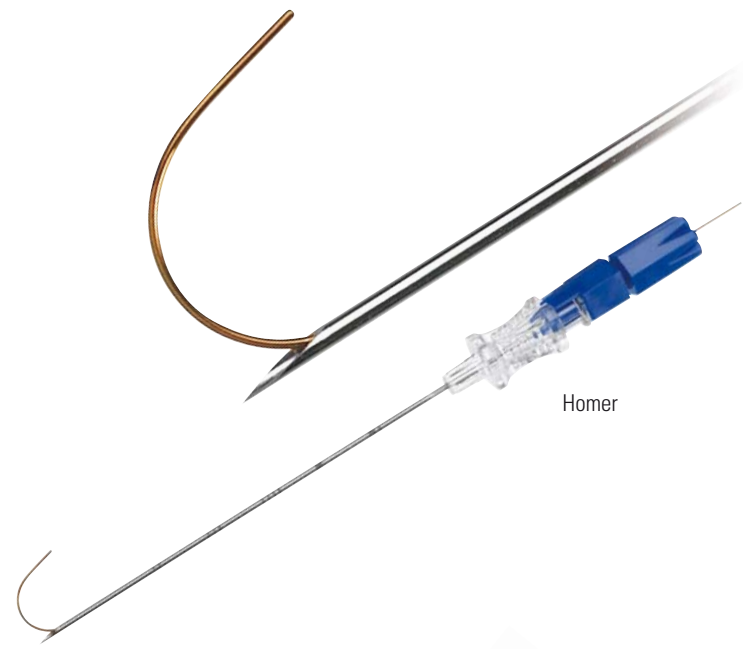

Homer

Homer Breast Localization Needle (Box of 10)

CATALOG NO.	DESCRIPTION
231030G	20G x 3 cm
231050G	20G x 5 cm
231075G	20G x 7.5 cm
231100G	20G x 10 cm
231125G	20G x 12.5 cm

Hawkins™ II is repositionable to achieve accurate placement.

- Available in two wire choices: FlexStrand™, a 19-strand cable, for maximum flexibility with minimal risk of transection, or a traditional hardwire, which minimizes flex during placement
- Both wires offer a side barb for temporary fixation during a confirming mammogram. The needle is easily repositioned, if necessary
- Easy-to-read centimeter markings facilitate accurate depth placement.
- Throated needle hub for easy wire loading
- Skin retention clip prevents forward migration of wire during patient transport
- Wire is marked in three places: I indicates wire at distal tip of needle; II indicates side barb deployment; and III indicates distal barb deployment
- Standard marking cannula on wire, beginning 2.5 cm from distal tip, serves as a landmark for the surgeon



Hawkins™ II

Hawkins™ I Breast Localization Needle (Box of 10)

CATALOG NO.	DESCRIPTION
251050	20G x 5 cm
251075	20G x 7.5 cm
251100	20G x 10 cm
251125	20G x 12.5 cm

Hawkins™ II Breast Localization Needle w/Hardwire (Box of 10)

CATALOG NO.	DESCRIPTION
242050	20G x 5 cm
242075	20G x 7.5 cm
242100	20G x 10 cm
242125	20G x 12.5 cm

Hawkins™ II Breast Localization Needle w/FlexStrand™ (Box of 10)

CATALOG NO.	DESCRIPTION
252050	20G x 5 cm
252075	20G x 7.5 cm
252100	20G x 10 cm
252125	20G x 12.5 cm