


Atrieve™ Vascular Snare Kit

Designed to fit a range of vessel sizes, the presence of three loops increases the probability of foreign body capture and provides a distinct advantage when working in confined spaces.

- Multiple, independent, short profile loops with tips on each loop allow for self-sizing to adjust to changing vessel diameter and complete vessel coverage
- Platinum strands provide excellent fluoroscopic visualization
- Shaft and loops made with nitinol wire for 1:1 torque control and kink-resistance
- 15° angled tip (6F and 7F sizes only) on catheter assists with navigation through tortuous anatomy
- Radiopaque band on catheter for increased visibility under fluoroscopy
- Kit includes an Atrieve™ Vascular Snare, delivery catheter, introducer and torque device

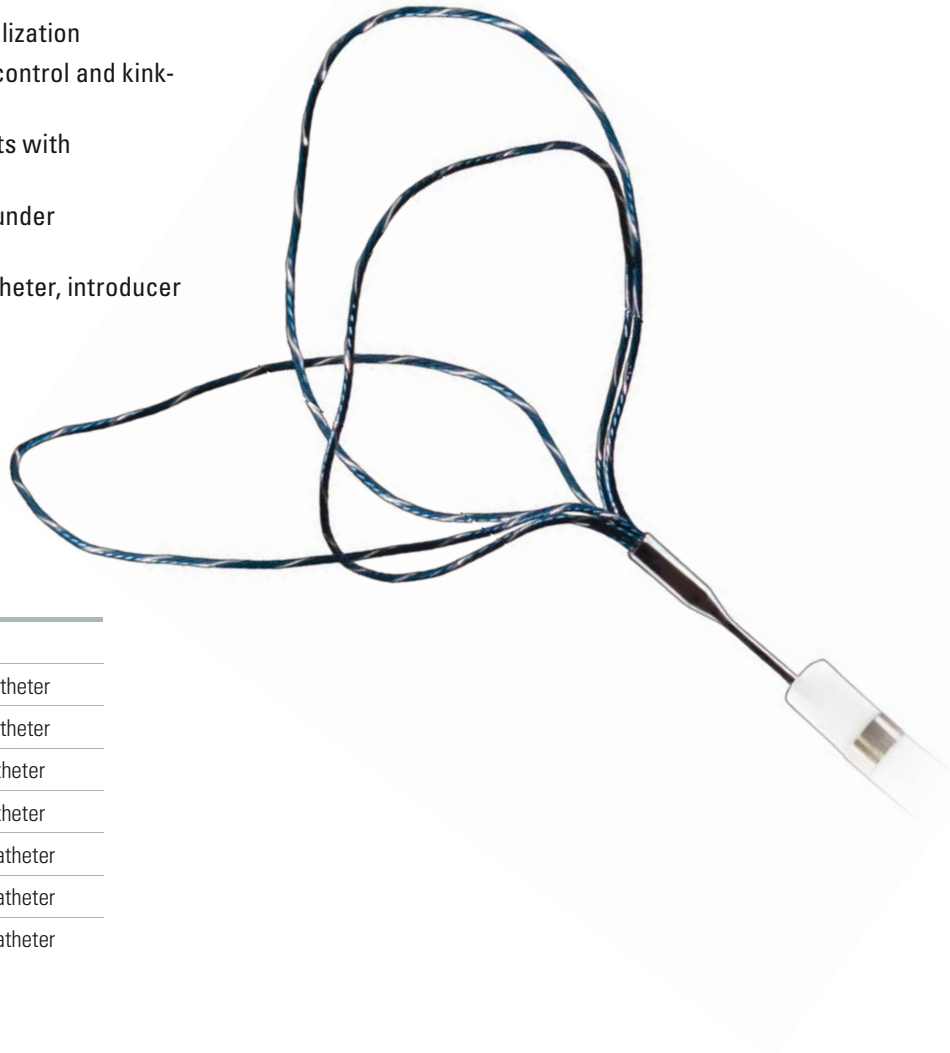

Atrieve™ Vascular Snare Kit (Each)

CATALOG NO.	DESCRIPTION
381003004	2-4 mm diameter x 175 cm snare, 3.2F x 150 cm catheter
381003008	4-8 mm diameter x 175 cm snare, 3.2F x 150 cm catheter
382006010	6-10 mm diameter x 120 cm snare, 6F x 100 cm catheter
382006015	9-15 mm diameter x 120 cm snare, 6F x 100 cm catheter
382006020	12-20 mm diameter x 120 cm snare, 6F x 100 cm catheter
382007030	18-30 mm diameter x 120 cm snare, 7F x 100 cm catheter
382007045	27-45 mm diameter x 120 cm snare, 7F x 100 cm catheter

