

New device for fixation of all drainage catheters


Reduce risk of catheter dislocation - without suturing

SKATER® FIX is a latex free catheter fixation device, consisting of a central thin carrier which has hydrocolloid adhesive and is covered with non woven medical grade adhesive. The soft foam support block and velcro prevent the catheter from kinking.

The catheter will be mechanically fixed with velcro, and minimal tape is involved directly on the catheter.

The velcro cover can be opened for puncture site inspection. In this way you have the possibility to look for infection or leakage whenever required.

- Latex free.
- Safe anchoring of drainage catheter up to one week.
- Allows inspection of puncture site without change of dressing.
- Easy to change.
- Minimal risk of catheter dislocation during change of fixation disk.
- Foam bridge prevents catheter kinking.
- Clean catheter with no tape residuals.
- Adjustable to catheters from 5 to 16 French.
- Hydrocolloid carrier will absorb body fluids from puncture site.
- Minimum risk of skin irritation.
- Disc follows patient contour and movements.

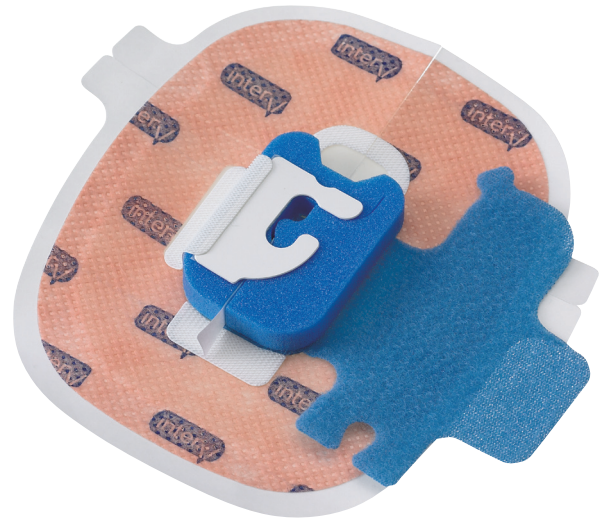

Special website with user guide and frequently asked questions: www.skaterfix.com

SKATER® FIX catheter fixation device

Box of 25

8025 00 001 Latex free catheter fixation device