

Non-Repositionable Breast Localization Needles

Hawkins™ III

- Available in two wire choices: FlexStrand™, a 19-strand cable providing maximum flexibility with minimal risk of transection, or a traditional hardwire that minimizes flex during placement
- Centimeter markings facilitate accurate depth placement
- Throated needle hub for easy wire loading
- Skin retention clip prevents forward migration of wire during patient transport
- Wire is marked in two places: I indicates wire at distal tip of needle; II indicates barb deployment
- Standard marking cannula on wire, beginning 2.5 cm from distal tip, serves as a landmark for the surgeon
- End hole design allows for dye injection or fluid aspiration
- Available with echogenic tip



Hawkins™ III

Hawkins™ III Breast Localization Needle (Box of 10)

CATALOG NO.	DESCRIPTION
244030	20G x 3 cm, hardwire w/echogenic tip
244050	20G x 5 cm, hardwire w/echogenic tip
244075	20G x 7.5 cm, hardwire w/echogenic tip
244100	20G x 10 cm, hardwire w/echogenic tip
244125	20G x 12.5 cm, hardwire w/echogenic tip

Hawkins™ III Breast Localization Needle w/FlexStrand™ (Box of 10)

CATALOG NO.	DESCRIPTION
253030	20G x 3 cm
253050	20G x 5 cm
253075	20G x 7.5 cm
253100	20G x 10 cm
253125	20G x 12.5 cm

Hawkins™ III Breast Localization Needle w/Hardwire (Box of 10)

CATALOG NO.	DESCRIPTION
243030	20G x 3 cm
243050	20G x 5 cm
243075	20G x 7.5 cm
243100	20G x 10 cm
243125	20G x 12.5 cm

For more information, please call the Angiotech Customer Care Center at 877 991-1110 (USA and Canada only) or 732 626-6466, option 1
 • Email interventional@angiotech.com • Visit www.AngioEduPRO.com, our educational website for healthcare professionals


Non-Repositionable Breast Localization Needles

Accura™ Breast Localization Needle (Kopans style) is a traditional, springhook hardwire with stiffener design.

- Skin retention clip prevents forward migration of wire during patient transport
- Wire is marked in two places: I indicates wire at distal tip of needle; II indicates barb deployment
- End hole design allows for dye injection or aspiration of fluid

Accura™ Breast Localization Needle (Kopans™ style) (Box of 10)

CATALOG NO.	DESCRIPTION
BLN2003	20G x 3 cm
BLN2005	20G x 5 cm
BLN20075	20G x 7.5 cm
BLN2010	20G x 10 cm
BLN2015	20G x 15 cm
BLN2103	21G x 3 cm
BLN2105	21G x 5 cm
BLN21075	21G x 7.5 cm
BLN2110	21G x 10 cm

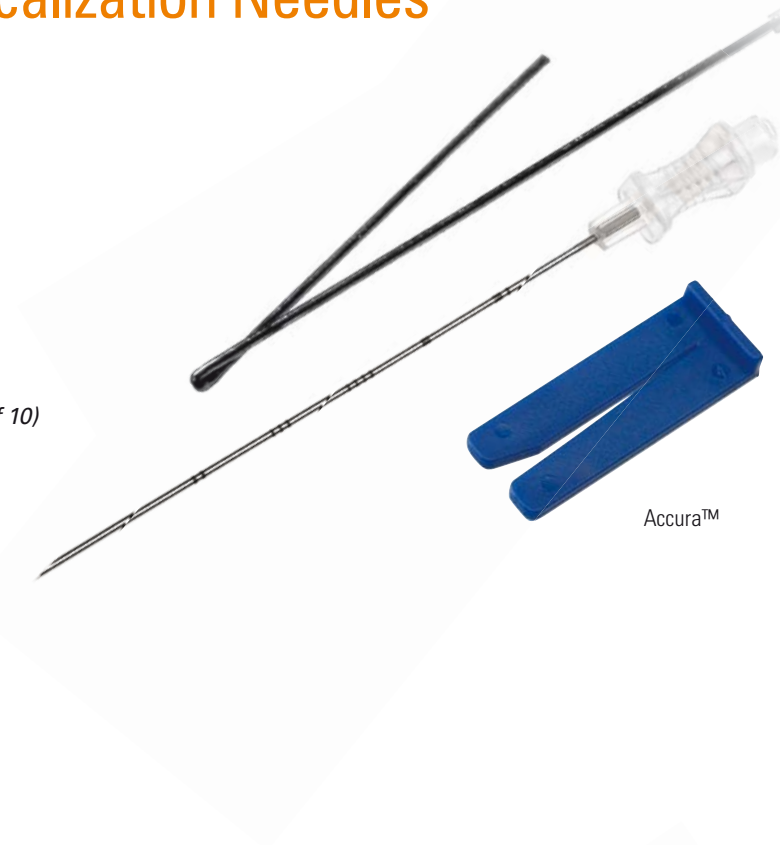


Accura™

Accura™ II Breast Localization Needle (Sadowsky™ Style) comes with a sterile separately packaged 23G stiffening cannula. Easily placed over the wire by the surgeon, the cannula indicates the course of the wire, provides a palpable landmark, and prevents accidental wire transection.

- Centimeter markings facilitate accurate depth placement
- Throated needle hub for easy wire loading
- Wire is marked in two places: I indicates wire at distal tip of needle; II indicates barb deployment

Accura™ II Breast Localization Needle (Sadowsky™ style) (Box of 10)

CATALOG NO.	DESCRIPTION
266050	20G x 5 cm w/23G stiffening cannula
266075	20G x 7.5 cm w/23G stiffening cannula
266100	20G x 10 cm w/23G stiffening cannula
266125	20G x 12.5 cm w/23G stiffening cannula

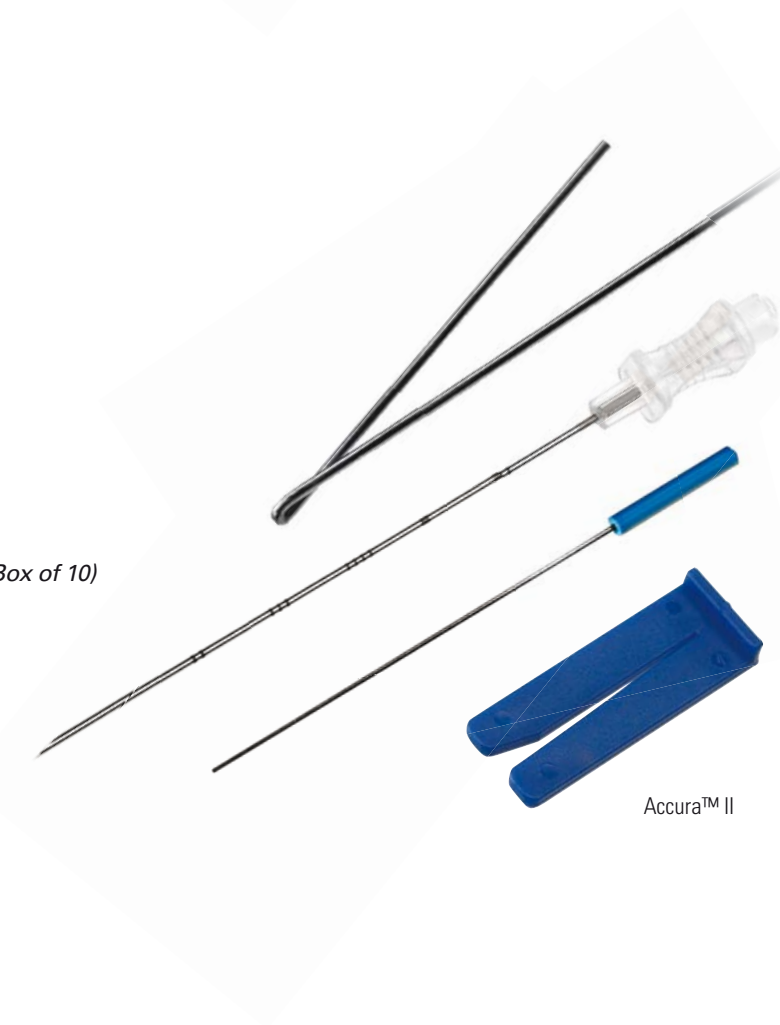

Accura™ II

Non-Positionable Breast Localization Needles

D Wire Breast Localization Needle

- A D-shaped wire which provides an increased cross sectional area aimed at reducing the risk of surgical transection. The tip of the hookwire is welded to increase anchoring strength
- Centimeter markings facilitate accurate depth placement. A throated needle hub gives easy wire loading
- End hole needle allows for dye injection or aspiration of fluid

D Wire Breast Localization Needle *(Box of 10)*

CATALOG NO.	DESCRIPTION
BLND2003	20G x 3 cm
BLND2005	20G x 5 cm
BLND20075	20G x 7.5 cm
BLND2010	20G x 10 cm
BLND2015	20G x 15 cm



D Wire