



# ROI-T™

TRANSFORAMINAL APPROACH IMPLANT



# ROI-T™

## TRANSFORAMINAL APPROACH IMPLANT



### Implant Features



- Minimal implant width and tapered tip facilitate insertion through a small opening.
- Chevron shaped teeth give all direction stability.
- Anatomic shape conforms to the anterior vertebral body.
- Open interior cavity allows for substantial graft amounts.
- Anterior and posterior side holes enable graft vascularization.
- Two lordosis options (5° and 9°) restore and maintain sagittal balance.
- Multiple heights (8mm to 16mm in 1mm increments) adapt to patient anatomy for a more secure construct.
- PEEK-Optima® material for biocompatibility and radiolucency.
- Tantalum markers for verification of implant positioning.
- Sterile packaging assures product quality and absolute traceability.



### Instrumentation

- Instrumentation enables easy access, insertion and placement.
- Shape-matched trials permit accurate fit assessment.
- Inserter fully supports the implant, providing a strong connection and placement maneuverability.
- Inserter's strong attachment also enables a superior assessment of trial stability.

### Dimensions

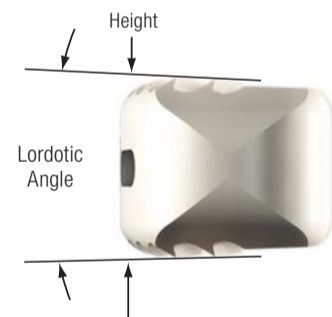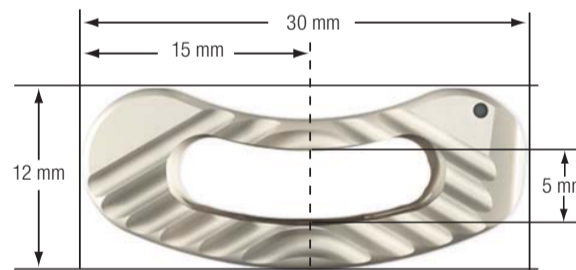


**Footprint (Width x Length):** 12mm X 30mm  
**Lordosis Options:** 5° and 9°  
**Height Options:** 8, 9, 10, 11, 12, 13, 14, 15, and 16mm

### Tantalum Markers



AP View

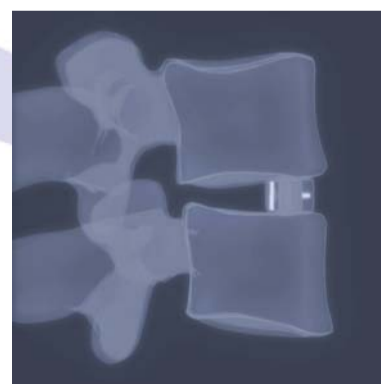

Lateral View





a passion for innovation

[www.ldrmedical.com](http://www.ldrmedical.com)

**France**

Technopôle de l'Aube BP 2  
10902 Troyes Cedex 9  
France  
+33 (0)3 25 82 32 63

**China**

Beijing Global Trade Center #36  
North Third Ring Road East,  
Unit 06, Level 19, Building A,  
Dongcheng District,  
Beijing, China, 100013  
+86 10 58256655

**Brazil**

Av. Pereira Barreto, 1395  
Torre sul - CJ 193 - Bairro Paraiso  
Santo André - São Paulo  
CEP: 09190-610  
Brazil  
+55 11 43327755

**United States**

4030 West Braker Lane, Suite 360  
Austin, Texas 78759  
512.344.3333

LDR, LDR Spine, LDR Médical, BF+, BF+(ph), Easyspine, Laminotome, MC+, Mobi, Mobi-C, Mobi-L, Mobidisc, ROI, ROI-A, ROI-MC+ and ROI-T are trademarks or registered trademarks of LDR Holding Corporation or its affiliates in France, the United States and other countries.